


JATI AUTOMASI ENTERPRISE

(Co. no 001773550-P)

No. 9B Jalan SG ¾, Taman Sri Gombak, 68100 Batu Caves, Selangor Darul Ehsan.
Tel : 03-61859326 Fax : 03-61871294

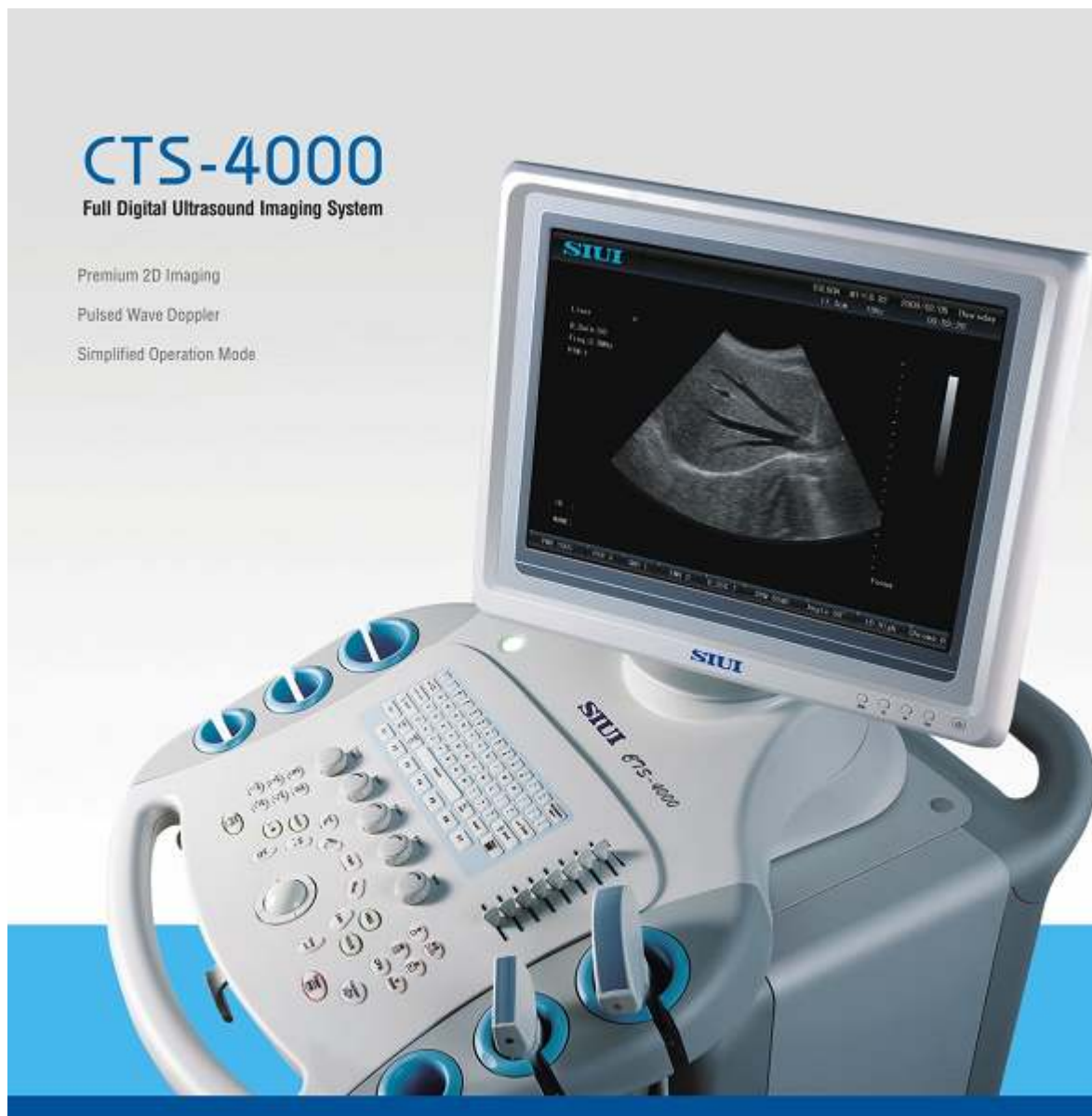

JATI AUTOMASI ENTERPRISE

(Co. no 001773550-P)

No. 9B Jalan SG ¾, Taman Sri Gombak, 68100 Batu Caves, Selangor Darul Ehsan.
Tel : 03-61859326 Fax : 03-61871294



Outstanding 2D imaging

The CTS-4000 is in itself a high-end trolley ultrasound system, and it is featured with highly clear image quality and strong functions.

High-Precision Digital Imaging Technology

- HiFi Digital Beam Forming (HDBF)
- Realtime Continuous Dynamic Focusing (RCDF)
- Realtime Dynamic Frequency Scanning (RDFS)
- Realtime Dynamic Aperture (RDA)
- Dynamic Realtime Apodization (DRA)
- High-Density Beamforming Scanning (HDBS)



Accurate beam-forming and signal processing, digital image acquisition and processing ensure images with clear-cut edge and no distortion.

Excellent Image Processing Technology

- Tissue Characterization Imaging
Based on characteristics of different tissues, the optimal parameters are preset in the system for acquiring ideal images easily.
- High Frequency Sampling
To deliver images with real representation of actual organs.



JATI AUTOMASI ENTERPRISE

(Co. no 001773550-P)

No. 9B Jalan SG ¾, Taman Sri Gombak, 68100 Batu Caves, Selangor Darul Ehsan.
Tel : 03-61859326 Fax : 03-61871294

Powerful Function

Five-frequency Probe

Each probe has five frequencies for selection. Finer images can be obtained according to different diagnosis needs.



Convertible Image Orientation


Images can be inverted not only up/down/left/right, but 90° continuously. The doctors may observe images as they wish, clear at a glance in clinical exchange.



Original (Liver) 180° (Uterus) 90° CW (Kidney) 90° CCW (Eyeball)

Unique High-Definition Zoom Function

When scanning images, the user can double acquire data from the region of interest, making the image finer and tissue resolution higher.

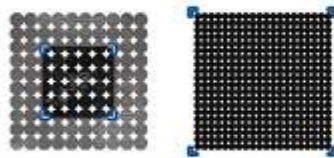

X1.0, sampling point: 100/cm² X2.0, sampling point: 400/cm²

Smart Tracing Function

When measuring area or circumference, if the trace is not good, press Backspace to return to the proper place and then continue tracing.



Complete M-Mode Function

- Editable screen layout for B-mode and M-mode:
Up/down or left/right display for selection.
- 2560-second super large capacity for M-mode cine.

Auto Freeze Function

To extend the service life of probes effectively.

User-defined Function Key

The user can assign frequently used functions on certain keys. For example, assign the BPD measurement function on F1 key. Press F1 and then you can access the function, which brings easy and convenient.

PIP (Picture-in-Picture) Function

The user can zoom in any region of interest, while maintaining the original diagnostic image for reference, to observe the zoom-in part clearly, thus diagnostic efficiency is improved.


IP One-key Optimization

One key for eight-parameter adjustment, easy image optimization.

Screen Saver Function

To protect the special ultrasound monitor and extend its service life.

Pulsed Wave Doppler (Option)

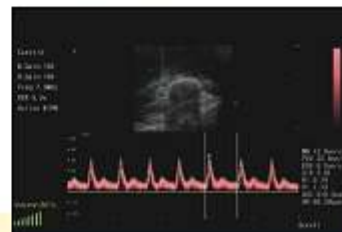

Measure artery flow in the palm with auto envelope, flow velocity S/D ratio RI PI acceleration HR etc.